


# paradigm™

A new paradigm  
of performance  
in endodontics

Paradigm BX endodontic  
files are treated with  
**LaserFlex technology.**

We use a **patented laser process** to selectively program the NiTi files with austenite (B2) and martensite (B19) phases across the section of the file. The resulting hybrid mixture of micro-structures in the files provides greater torsional strength, durability and cutting edge retention with **20% more flexibility** compared to any leading heat-treated Niti files.

Some heat-treated files provide flexibility by tempering the microstructure of the Nitinol base material - but this flexibility gives up sharpness and hardness. With the **LaserFlex** technology we use for the **Paradigm BX files**, we have achieved a highly flexible material that will stay sharp and last longer.

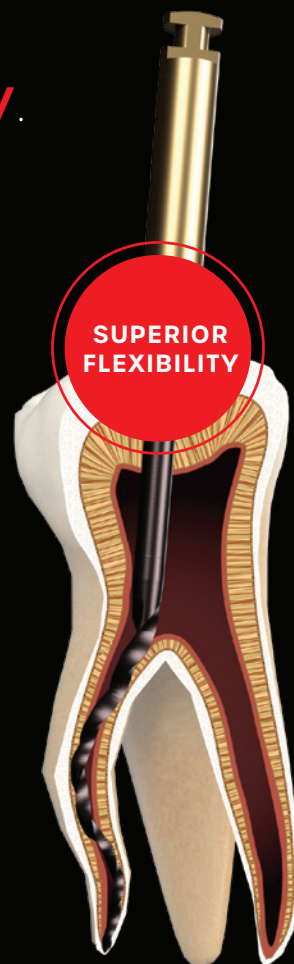

#### CONFIDENCE

Paradigm BX files are more resistant to separation from torsional overload than the leading heat-treated multi-taper file.



#### EFFICIENCY

Paradigm BX files have a hard cutting edge to stay sharp longer for cutting efficiency and superior cleaning.



#### FLEXIBILITY

Paradigm BX files are 20% more flexible than other leading heat-treated multi-taper file for improved performance in difficult canals.



#### ENDURANCE

Paradigm BX files are laser processed to a hybrid microstructure for torsional resistance without giving up durability in cyclic bend fatigue life.

# paradigm BX™


LIMITED TIME  
INTRODUCTORY OFFER

**50% OFF**

UNTIL AUGUST 31, 2021

## FEATURES

- Patented LaserFlex process
- Hybrid BX microstructure
- Soft martensitic flexible core
- Hard austenitic cutting edge
- Multiple-taper design
- Non cutting tip
- Fits all standard handpieces

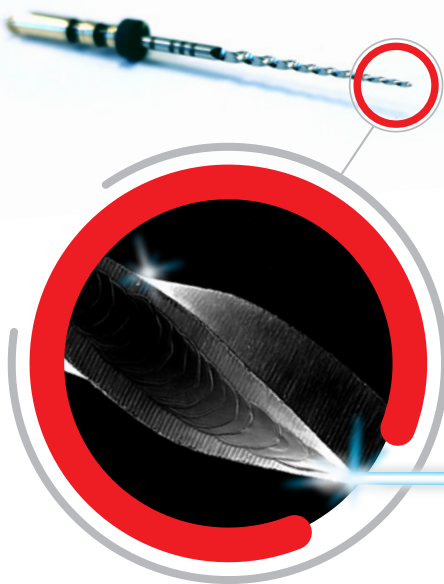

SIZE	LENGTH	TIP DIA.	TAPER
AX	19 mm	018	Multi-taper
A1	25 mm	017	Multi-taper
A2	25 mm	018	Multi-taper
M1	25 mm	020	Multi-taper
M2	25 mm	025	Multi-taper
M3	25 mm	030	Multi-taper
M4	25 mm	040	Multi-taper
M5	25 mm	045	Multi-taper

The **LaserFlex** process makes **Paradigm BX** the most technologically advanced file ever.

With more flexibility and longer-lasting surfaces, these files **increase confidence, comfort, and patient care** with no interruption to your workflow.

Contact us today to purchase the advanced **Paradigm BX** NiTi files and take advantage of our **introductory offer!**

☎ 1.833.256.0404   ✉ sales@paradigmendo.com   ▶ paradigmendo.com